



SampleManager LIMS[™] OData Service

Thermo Scientific[™] SampleManager[™] LIMS 21.1 OData service allows external applications or web publishing services to obtain data from your SampleManager LIMS[©] instance.

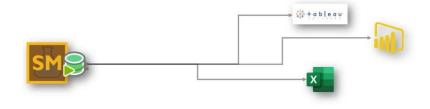

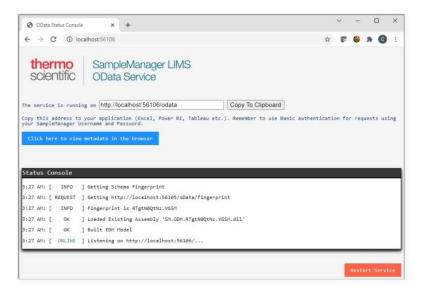

It differs from direct database access as it utilizes the enhanced data governance of SampleManager LIMS[©], honors roles, limits which fields can be seen, and provides custom functions. External applications can only see the query definitions provided by the service.

SampleManager LIMS[©] OData Service utilizes the OData queries to define the selection criteria on rows from a SampleManager entity type, the columns of the selected data that are available through the SampleManager LIMS[©] OData service, and an optional set of functions to process the selected data. The Service Uses Its Own Dedicated Windows / Console Service To Safely Deliver Specific Data To External Applications.

Being more flexible than other REST-based web services, it provides a uniform way to describe the data and the data model for easy interoperability between data sources, applications, services, and clients. The SampleManager LIMS[®] OData Service is supported by applications such as Excel, Power BI, Tableau, and many more.

ODATA QUERIES

>	CurrentSamples	
Description	Current samples (set 7 days) with results	
Entity Type	SHIP_TEST_RESULT	
Group		
Modified On	02/08/2022 10:27	
Modified By	Listen Nereier	
Modifiable	2	Published

Used by permission of Thermo Fisher Scientific, the copyright owner